
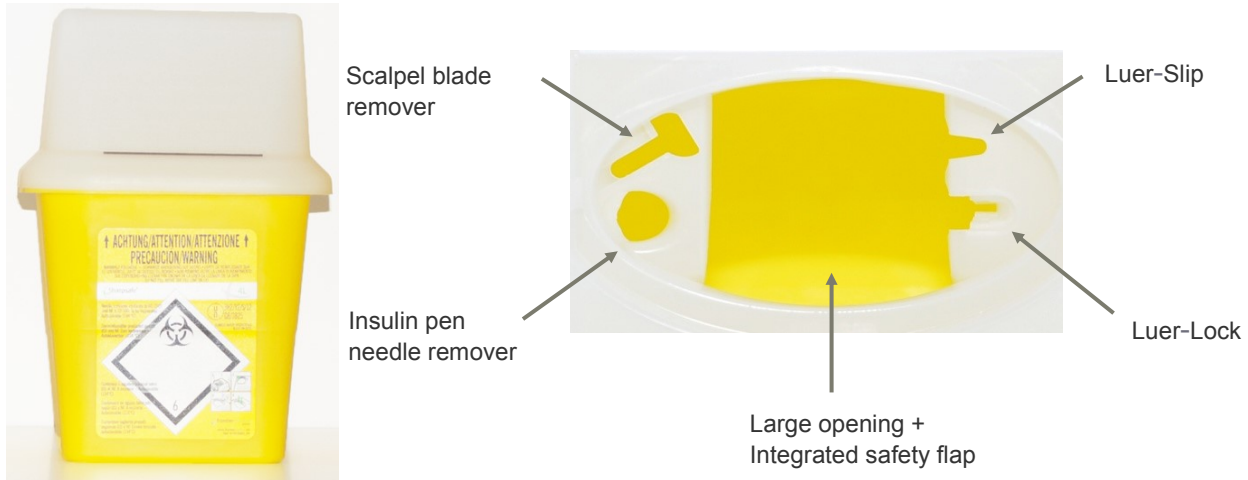


Sharpsafe[®] 4 L Medical Waste - The classic

SHA-41055720

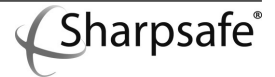

Features

- The container consists of two parts : a yellow body and translucent lid
- Safe, quick and simple 4 pressures assembly with 4 audible « click »
- Disconnection for all types of needle even for the scalpel blade and insulin pen
- 2 closure positions on the flip top lid: temporary and final closure
- Reopening of the temporary closure with a simple finger pressure on the lid
- Possibility to block to lid in the open position
- Safety flap for extra protection
- Clearly visible fill line on the translucent lid
- 6cm safety area between the fill line and the opening of the container
- In-mould label which cannot be removed, even by autoclaving
- Handle for a safe transport
- Can be used in a "single hand" way
- Accessories: Trays, stands, hooks and brackets available
- Instruction of use on the label: German, French, Italian, Spanish and English
- Autoclaving valid to 134°C during 18 minutes for an empty and open container
- Conform to the new ISO Norme 23907:2013, TRBA 250, NFX30-500 and BS7320 in terms of resistance to penetration, drop impact and liquid resistance. ADR and NF302 certificated. Holder of the brand NF302 (details see reverse)
- Origin of the product : European Union

HUBERLAB. AG
Industriestrasse 123
4147 Aesch

T 061 717 99 77
F 061 711 93 42

www.huberlab.ch
info@huberlab.ch








Technical information

- Material 100 % Polypropylen (Without Latex)
- Max. nominal volume 4 L
- Max. filling capacity 3 L
- Dimensions opening 55 x 50 mm
- Dimensions flip top closed 245 x 175 x 175 mm (H x L x W)
- Weight 0,320 kg

Packaging

- Quantity per box 50
- Dimensions box 416 x 840 x 355 mm (H x L x W)
- Weight 16,8 kg
- Volume 0,125 m³
- Box per pallet 14
- Pallet dimensions 226 x 120 x 80 cm
- Palette weight 255 kg

Conformity

	Die ISO Norm 23907:2013 describes the requirements and the test methods concerning the production of the sharps containers with the objective to protect the personal against sharps injury. The ISO norm is through Kitemark certified. Sharpsafe was the first brand to obtain this certification.
	The NF X 30-500 is the latest Sharps-guideline in Europe and the most rigorous European standard in terms of resistance to penetration, drop impact, perforation and liquid resistance
	The NF302 (Norme française) ensures the continued existence of the standards defined by NF X 30-500. The 'laboratoire National d'Essais' (French institute for Standards) checks regularly the quality of the products on the market
	In conformity with ADR and UN3291: International norm for the transport of dangerous material and infectious waste
	The TRBA 250 is a German guideline for the treatment of biological, dangerous and infectious waste
	Quality management system : ISO 9001:2008 et ISO 13485:2003 Environmental standard: ISO 14001:2004
	Autoclave at 134°C during 18 minutes